

Glassware Washer Buying Guide


Undercounter & Freestanding/Mobile Units



436601



UCW-HA10101

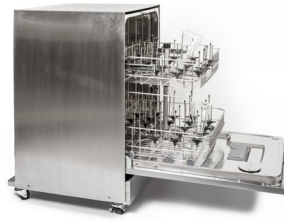

Start Here

Do you need a single unit (6ft³ volume with 2 racks) or are you interested in two, vertically stacked units (12ft³ volume with 4 racks)

Two Vertically Stacked Units
See Vertical SpaceSaver™ units on reverse side

Single Unit
Do you need the unit to fit under a counter?

Undercounter
Do you need the washer to hold narrow neck glassware such as volumetric or Erlenmeyer flasks?

Freestanding/Mobile
Do you need the washer to hold narrow neck glassware such as volumetric or Erlenmeyer flasks?

No
We recommend **UCW-HA00101**, Standard Freestanding & Mobile Glassware Washer

Yes
You need a Spindle Rack Ready* washer.
Is your incoming DI water pressure 18 psi or greater?

No
Do you need a DI Rinse Function?

Yes
You need a Spindle Rack Ready washer.
All units include a DI rinse function. Is your incoming DI water pressure 18 psi or greater?

No
Do you need a DI Rinse Function?

Yes
Is the incoming DI water pressure 18 psi or greater?

Yes
You need **UCW-HA20001**, Undercounter Glassware Washer with DI Rinse, Spindle Rack Ready*

No
For incoming DI water pressure under 18 psi, you need a DI booster pump. Try **UCW-HA20001**, Undercounter Glassware Washer with DI Rinse, Spindle Rack Ready* and **446508**, Booster Pump

Yes
Is the incoming DI water pressure 18 psi or greater?

No
We recommend **UCW-HA00001**, Standard Undercounter Glassware Washer

Yes
We recommend **UCW-HA20101**, Freestanding & mobile Glassware Washer with DI Rinse and Spindle Rack Ready*

No
For incoming DI water pressure under 18 psi, you need a DI booster pump. Try **UCW-HA20101**, Freestanding & Mobile Glassware Washer with DI Rinse, Spindle Rack Ready* and **446508**, Booster Pump

Yes
We recommend **UCW-HA10101**, Freestanding & Mobile Glassware Washer with DI Rinse

No
We recommend **UCW-HA10101**, Freestanding & Mobile Glassware Washer with DI Rinse and **446508**, Booster Pump

Yes
We recommend **UCW-HA10001**, Undercounter Glassware Washer with DI Rinse

No
For incoming DI water pressure under 18 psi, you need a DI booster pump. Try **UCW-HA10001**, Undercounter Glassware Washer with DI Rinse, Spindle Rack Ready* and **446508**, Booster Pump

*Spindle Racks Sold Separately:

100005409: 2 pc. Spindle Rack Set - full bottom rack, 1/2 open / 1/2 spindle top rack

436601: Spindle Rack, Bottom Only (For Spindle Rack Ready Washers Only)

4370802: 2 pc. Spindle Rack Combo Set (For Spindle Rack Ready Washers Only)



Glassware Washer Buying Guide

Vertical SpaceSaver™ Units



Start Here

Do you need a single unit (6ft³ volume with 2 racks) or are you interested in two, vertically stacked units (12ft³ volume with 4 racks)

Single Unit

See Undercounter & Freestanding/Mobile units on reverse side

Two Vertically Stacked Units

Vertical SpaceSaver™ units feature two independent washers securely stacked in a stainless steel housing. Do you need the washer to hold narrow neck glassware such as volumetric or Erlenmeyer flasks?

Yes

You need a Spindle Rack Ready* washer. Is your incoming DI water pressure 18 psi or greater?

No

Do you need a DI Rinse Function?

Yes

Do you need the same Spindle Rack Ready* features for top unit?

No

For incoming DI water pressure under 18 psi, you need a DI booster pump. We recommend also purchasing **446508**, Booster Pump

Yes

All bottom units with DI rinse are also Spindle Rack Ready,* but may be used with included standard racks. Is the incoming DI water pressure 18 psi or greater?

No

If bottom unit is a standard water washer, top unit must be the same. We recommend **VSW-HA0001**, Vertical SpaceSaver™ Glassware Washer for Use with Tap Water

Yes

We recommend **VSW-HA22001**, Vertical SpaceSaver™ Glassware Washer with DI Rinse, Top and Bottom Spindle Rack Ready*

No

We recommend **VSW-HA12001**, Vertical SpaceSaver™ Glassware Washer, Top Standard Rack, Bottom Spindle Rack Ready*

Yes

When bottom unit has DI Rinse Function, top unit must have the same DI Rinse Function, but does not have to be Spindle Rack Ready.* We recommend **VSW-HA11001**, Vertical SpaceSaver™ Glassware Washer with DI Rinse, Top and Bottom Standard Racks

No

For incoming DI water pressure under 18 psi, you need a DI booster pump. We recommend **VSW-HA11001**, Vertical SpaceSaver™ Glassware Washer with DI Rinse, Top and Bottom Standard Racks and **446508**, Booster Pump

*Spindle Racks Sold Separately:

100005409: 2 pc. Spindle Rack Set - full bottom rack, 1/2 open / 1/2 spindle top rack

436601: Spindle Rack, Bottom Only (For Spindle Rack Ready Washers Only)

4370802: 2 pc. Spindle Rack Combo Set (For Spindle Rack Ready Washers Only)

Contact SP Hotpack:

SP Corporate Office
935 Mearns Road
Warminster, PA 18974 USA
U.S.: +1-800-523-2327
Intl: +1-215-672-7800
sales@spscientific.com

www.sp-hotpack.com

SP-04076-RevA

